
**Oral-B**

**iO™**

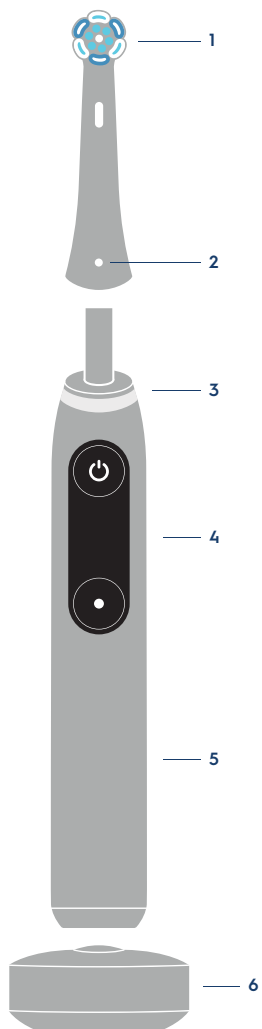9

10

11

12

13

14

15

16

17

18

18

19

20

21

22

23

24

25

26








27

28

29

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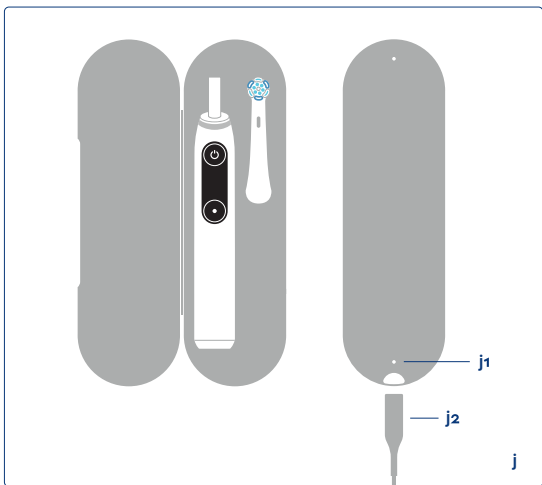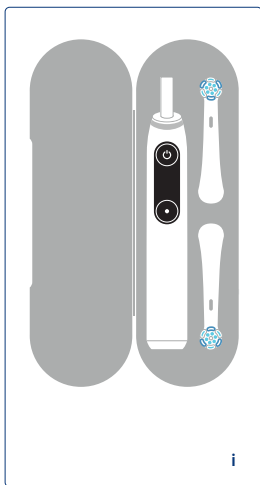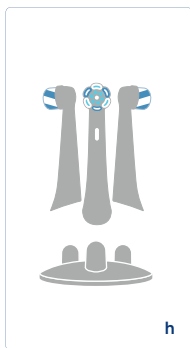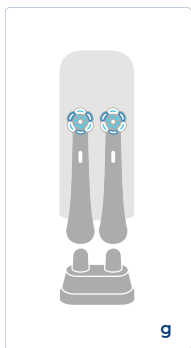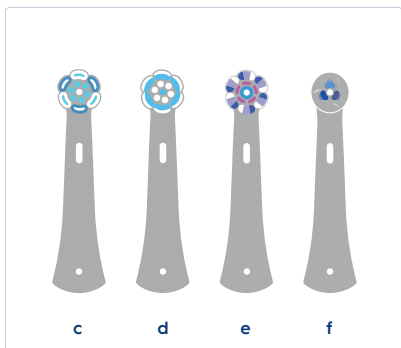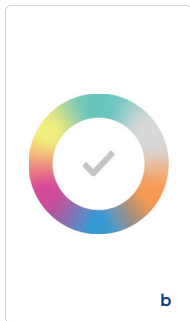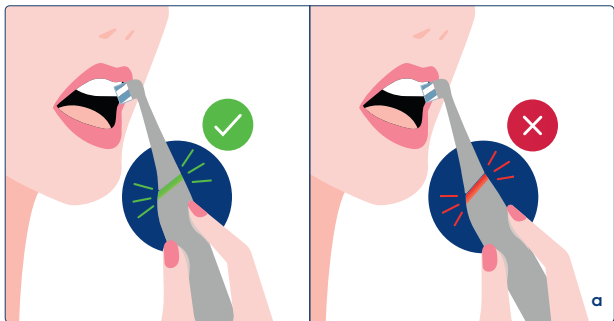## Accessories (depending on model) and features:

- a) Pressure Sensor features
- b) Customise colour of Light Ring via Oral-B App
- c) Brush head «Ultimate Clean»
- d) Brush head «Gentle Care»
- e) Brush head «Radiant White»
- f) Brush head «Specialised Clean»
- g) Brush head holder with cap
- h) Brush head holder stand
- i) Premium travel case
- j) Charging travel case
  - 1 Charging light
  - 2 Power supply (Smart Plug)
- k) Download Oral-B App and connect to *Bluetooth*<sup>®</sup> (see back page)

**Note:** Content may vary based on the model purchased.

# English

## Rechargeable toothbrush.

(UK)(IE)

### Welcome to Oral-B!

Before operating your Oral-B iO toothbrush, please read these instructions and save this manual for future reference.

### IMPORTANT

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- If the product/cord is damaged, take it to an Oral-B Service Centre.
- Only use charger provided with your appliance.
- Recommended ambient temperature for charging is 5°C to 35°C.  
Do not expose the appliance to temperatures higher than 50°C for extended periods of time.

### WARNING

- Do not place the charger, smart plug or charging travel case in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- This appliance contains a battery that can be replaced. Please note that battery replacement is only available in certain countries. Before considering battery replacement, visit [www.service.oralb.com](http://www.service.oralb.com) for information and safety precautions.
- Usage by children under age 3 is not recommended.
- Do not chew or bite on brush. Do not use as a toy or without adult supervision.
- Do not chew or bite on nozzle. Do not use as a toy or without adult supervision.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- For hygienic reasons we recommend the use of individual brush handles and brush heads for each family member.
- Small parts may come off, keep out of reach of children.
- Clean brush head properly after each use (see «Cleaning guidelines» paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.

To avoid brush head breakage which may generate a small part choking hazard or damage teeth:

- Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fit properly. Never use without a brush head.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Replace brush (head) every 3 months or sooner if brush head becomes worn.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- For tooth brushing only.
- Periodically check the entire product/cord/accessories for damage. A damaged or non-functioning unit should no longer be used.
- Do not modify or repair the product. This may cause fire, electric shock or injury.
- Never insert any object into any opening of the appliance.
- Avoid sharp bending, jamming or squeezing of the cord.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- Only use the charger and/or special cord set (power supply) with integrated Safety Extra Low Voltage power supply provided with your appliance. Do not exchange or tamper with any part of it, otherwise there is risk of an electric shock. If charging stand/charging travel case is marked with  $\text{D-292-C 492-XXXX}$ , only use Braun/Oral-B special cord set with this identification.
- This toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

### Important Information

- This toothbrush, charger and pouch (accessory) contain NdFeB magnets and generate magnetic fields. If you use a medical device or have an implanted medical device, consult with the medical device manufacturer or a physician, if you have questions about safely using your Oral-B iO toothbrush.
- To avoid electromagnetic interference and/or compatibility conflicts, disable *Bluetooth®* on your toothbrush handle before using it in restricted areas.

- Your Oral-B toothbrush can be used with your smartphone (find details under «Connecting your toothbrush to your Smartphone»).
- Disable *Bluetooth* by entering the «Settings» menu on the interactive display (4) → *Bluetooth* → disable. Follow the same procedure to enable *Bluetooth* again.
- Persons with pacemakers should always keep the toothbrush more than 15 centimeters (six inches) from the pacemaker while turned on. Anytime you suspect interference is taking place, disable *Bluetooth*.

### Specifications

For voltage specifications refer to the bottom of the charging unit.

Noise level: 64 dB (A)

## First Use

1. Please wake up your iO toothbrush by pressing the on/off button (7) and follow the instructions on the interactive display (4) to get started for first time setup.
2. If not further instructed on the display, you may toggle with the mode button (8) and select with the on/off button (7).
3. You may change the settings via the «Settings» menu (18) at any time.

## Charging and Operating

- To charge the toothbrush battery place the toothbrush on the plugged-in charger (6).
- Once the toothbrush is placed on the charger, the charge level is shortly indicated on the interactive display (4). The light ring (3) fades in and out during charging. Once it is fully charged the light ring turns off. A full charge typically takes 3 hours.
- You may check the battery capacity at various time points:
  - Before brushing by pressing the mode button (8) for 2 seconds.
  - After brushing within the brushing feedback screen (28).
  - When display turns off.
- If the battery is running low, the interactive display indicates the remaining battery capacity and the light ring (3) blinks in red. Once the battery is empty, the motor stops; it will need at least a 5-minute charge for one usage.
 

**NOTE:** In case the battery is empty (display off), charge for at least 30 minutes.

- You can always store the toothbrush on the plugged-in charger to maintain it at full charge; a battery overcharge is prevented by the device. However, for environmental purposes we recommend to keep charger in the power socket during charging only.

**Note:** Store the toothbrush at room temperature for optimal battery maintenance.

**CAUTION:** This toothbrush contains a Li-Ion battery. Do not expose the toothbrush to temperatures higher than 50 °C and lower than 5°C. In case the toothbrush is exposed to temperatures outside this range (e.g. during airplane travel) the brush will automatically switch off. You can re-activate your brush by putting it on the charger.

## Using your iO toothbrush

### Start brushing

1. Place the brush head on the toothbrush handle. A small gap should remain between handle and brush head.
2. Activate the display by picking up the toothbrush, lifting it off the plugged-in charger or pressing the on/off button (7) (display turns on).
3. Wet brush head and apply toothpaste. To avoid splashing, guide the brush head to your teeth before switching the toothbrush on by pressing the on/off button.
 

**Note:** Display turns off after some time. You may press the on/off button to activate display again.
4. Guide the brush head from tooth to tooth.
5. Brush the outsides, the insides and the chewing surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you.
6. After brushing (> 30 sec) the display provides feedback on your brushing behaviour. (25)

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

If you have sensitive teeth and/or gums, Oral-B recommends using the «Sensitive» or «Super Sensitive» modes (optional in combination with an Oral-B «Gentle Care» brush head (d)). Availability of modes may vary depending on model purchased.

### Toothbrush Timer (Stutter timer)

The iO toothbrush is equipped with timer features, which help you to achieve the professionally-recommended 2 minutes brushing time. The timer signals with a short stuttering sound at 30-second intervals to move to the next quadrant of your mouth. A long stuttering sound indicates the professionally-recommended 2 minutes brushing time.

The timer memorises the reached brushing time, even when the handle is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or placing the handle on the plugged-in charger.

### Brushing Reward

After completing a 2-minute brushing session, you will be rewarded with a starry-eyed smiley face (26). When brushing for 2 minutes and most of the time with ideal pressure, you will be rewarded with a starry-eyed smiley face and fireworks.

## Pressure Sensor / Toothbrush Light Ring

Your iO toothbrush is equipped with a light ring (3) that has various functions indicated via colours:

### Pressure Sensor (a)

- White = Low pressure indicator
- Green = recommended pressure indicator
- Red = high pressure indicator

If the Pressure Sensor illuminates in white, we recommend applying more pressure. The green light indicates recommended brushing pressure. If too much pressure is applied, it lights up in red reminding you to reduce pressure. During high pressure the brush head movement will be reduced (in «Smart Adapt», «Daily Clean», «Intense» and «Sensitive» modes).

Periodically check the operation of the Pressure Sensor by pressing moderately on the brush head during use.

**Note:** Pressure Sensor is deactivated during «Tongue Clean» mode (availability depending on model).

### Other Light Ring indicators

- Red = Low battery indicator
- Blue = *Bluetooth* pairing indicator
- Orange = Error indicator (follow display messages)

- Choose your favourite colour (20) via the «settings» menu.

## Brushing Feedback

After your brushing session (session min. 30 sec) your iO toothbrush provides feedback on your brushing behavior and toothbrush status. Toggle through the different screens by pressing the mode button (8). The availability of brushing feedback depends on model purchased:

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**Note:** You may customise the brushing feedback screens via the Oral-B App.

### Brushing modes (depending on model)

Your iO toothbrush offers different brushing modes indicated on the display:

«**Smart Adapt**» – Optimal cleaning and gum care through automatic intensity adaptation.

«**Daily Clean**» – Standard mode for daily cleaning

«**Sensitive**» – Gentle, yet thorough cleaning for sensitive areas

«**Whiten**» – Polishing for occasional or everyday use

«**Gum Care**» – Gentle massage of gums

«**Intense**» – For an extraordinary clean feeling (operated with a higher speed)

«**Super Sensitive**» – Super gentle cleaning for extra sensitive areas

«**Tongue Clean**» – Tongue cleaning for occasional or everyday use

**Note:** The sequence and availability of brushing modes may vary depending on model purchased.

When operating the «Tongue Clean» mode, we recommend using a «Gentle Care» brush head (d). You can brush your tongue either with or without toothpaste. Brush the whole area of your tongue systematically in gentle movements. The recommended brushing time is 30 seconds; the motor stops after 30 seconds.

#### Mode setting:

1. You may choose a mode by toggling through the mode menu via the mode button (8).
2. Select your preferred mode by pressing the on/off button (7). Your next brushing session will automatically start in your pre-selected mode.
3. You may also change modes during the brushing session.

**Note:** You may also customise mode settings via Oral-B App.

## Interactive Display

Activate the interactive display (4) by picking-up the toothbrush, lifting it off the charger or pressing the on/off button (7).

- (10-17) «Mode selection» – toggle through the mode menu by pressing the mode button and select your favourite mode via the on/off button.
- (18) «Settings» – toggle through all modes until you reach the «Settings» menu; to enter press the on/off button
- (19) «Brush Head replacement reminder» (indicates status of the brush head/ perform reset or enable / disable brush head replacement reminder)
- (20) «Light ring» - select/change the colour of the light ring
- (21) «Bluetooth» manage *Bluetooth*® settings (enable, disable)
- (22) «Language» - select/change the language
- (23) «Restart Brush» - restart toothbrush and revert settings to default
- (24) «Back» - return to mode menu

## Brush heads

Oral-B offers you a variety of different Oral-B iO brush heads (c/d/e/f) that only fit your Oral-B iO toothbrush handle. These unique brush heads are dentist-inspired to leave your mouth feeling clean and refreshed.

Brush head «**iO Ultimate Clean**» - for 360° Advanced Deep Cleaning

- Thousands of bristles angled at 16° for superior plaque removal
- Spiralled bristle pattern for better adaptability
- Dual-length bristles for cleaning even in hard-to-reach areas

Brush head «**iO Gentle Care**» - for Healthier Gums\* as of Day 1

- 4000+ soft bristles for effective cleaning
- Curved bristle pattern for gentle gliding effect
- Ultra-thin adaptable bristles – even for sensitive teeth & gums

\*vs. a regular manual toothbrush

Brush head «**iO Radiant White**» - for Whiter Teeth\* as of Day 1

- Thousands of densely arranged bristles for thorough cleaning
- Curved & angled outer bristles for better reach in between teeth
- 5 polishing petals to remove everyday stains (like the ones from coffee, wine, or smoking)

\*by removing surface stains

Brush head «**iO Interdental Clean**» - for Extra Clean in Hard-To-Reach-Areas

- Dual-length high-density bristles clean in between teeth
- Outer bristles angle in to support central bristles
- Use as Step 2 after your usual brushing routine – even implants & braces

Most iO Oral-B brush heads feature **SMART BRISTLES** to help you monitor brush head replacement need. The bristle colour will fade to white with usage, indicating the potential need to replace your brush head. If the bristles splay before the colour recedes, you may be exerting too much pressure on teeth and gums or scrubbing while brushing. In general, we recommend replacing the brush head every 3 months.

#### Brush head replacement reminder

The toothbrush is equipped with a «Brush Head replacement reminder» indicating the dentist-recommended replacement of the brush head after 3 months / 90 days. The time elapsed from the last

reset can always be checked in the brushing feedback (25) after the brushing session or in the «settings» menu. To reset the «Brush Head replacement reminder» follow the instructions on the display or perform the reset in the «settings» menu (19). You can also disable the screen via the «settings» menu.

**Note:** You may also customise the brush head replacement settings via the Oral-B App.

## Connecting your iO toothbrush to your Smartphone/ Oral-B App

The Oral-B App offers various features (depending on model purchased):

- **Personalise your toothbrush settings** (e.g. choose your preferred light ring colour, adjust stutter timer, brush mode sequence or brush head replacement reminder settings, enable/disable brushing feedback screens or select additional screens)
- **Track your brushing behaviour** (e.g. brushing time and applied pressure)
- **Real-time brushing recognition** to help you brush every area for the right amount of time.

The Oral-B App is available for smartphones running iOS or Android™ operating systems. Check [app.oralb.com](http://app.oralb.com) for compatible smartphones and country availability of the Oral-B App. It can be downloaded for free from App Store<sup>(SM)</sup>/Google Play store™ or you may scan the QR code on the back page of this manual (k).

1. Launch the Oral-B App. It will guide you through the entire *Bluetooth* pairing procedure. You may also find it under the App's settings (More → Your brush).  
**Note:** Oral-B App functionality is limited if *Bluetooth* is deactivated on your smartphone. (For directions please refer to your smartphone's user manual.)
2. First pairing is indicated via the blue light ring. Any App instructions will be displayed on your smartphone.
3. First pairing needs to be done via Oral-B App. The toothbrush is automatically connected with your smartphone for future brush events.
4. You may change *Bluetooth* settings on your toothbrush via the display → Settings → *Bluetooth* enable/disable.
5. Keep your smartphone in close proximity (within 5m distance) when using it with the


toothbrush. Make sure the smartphone is safely positioned in a dry place.

**Note:** Your smartphone must support *Bluetooth* 4.2 (or higher) / *Bluetooth* Smart in order to pair your toothbrush.


**CAUTION:** Follow your smartphone's user manual to ensure your phone/ charger is designed for use in a bathroom environment.

## Cleaning Guidelines

1. After each brushing, rinse your brush head under running water with the handle switched on.
2. Remove the brush head and clean the handle and brush head separately. Make sure water runs inside the top of the handle while the handle is switched on, as well as inside the brush head.
3. Wipe all parts dry and let them air-dry before reassembling your toothbrush again.

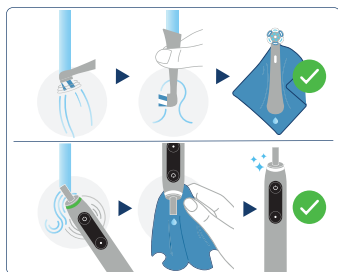4. Unplug the charger before cleaning.
5. Clean the charger regularly with a damp cloth.
6. Please check individual cleaning instructions for each accessory in its «Accessories» paragraph.
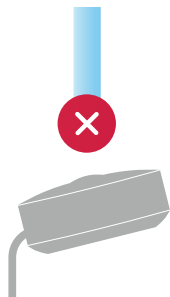
## Accessories (depending on model)

### Charging Travel Case (j)

- To charge your toothbrush (5) place it in the charging travel case (j) and connect it to an electrical outlet, using the extra power supply (j2) provided. The power supply may be used for all global mains voltage ranges (100 – 240V).
- The flashing light (j1) of the charging travel case indicates that the toothbrush is charging. Once the toothbrush is fully charged the light turns off. A full charge in the charging travel case typically takes 6 hours.
- The outer and inner case should be cleaned with a damp cloth only. Make sure that the toothbrush and brush head are completely dry before storing in the travel case.

**Note:** Clean stains on travel case surfaces immediately. Store travel case in a clean and dry place.

### Premium Travel Case (i)

Store the clean and dry iO toothbrush/ brush heads in the premium travel case (i). The outer and inner case should be cleaned with a damp cloth only.

**Note:** Clean stains on travel case surfaces immediately. Store travel case in a clean and dry place.

### Brush head holder stand (h) / Brush head holder with cap (g)

You may store your cleaned and dry brush heads on the pins of the brush head holders. Make sure your brush heads are dried before placing them in the holder and/or closing the cap. Brush head holders and cap are dishwasher-safe.

### Zipper Pouch / Magnetic Pouch

You may store your cleaned brush heads and toothbrush in the pouches. Make sure your brush heads and toothbrush are dried before placing them in the pouches. The outer and inner surfaces can be cleaned under running water. Store pouches in a clean and dry place.

**Note:** Clean stains on pouch surfaces immediately.

Subject to change without notice.

## Environmental Notice

Product contains a battery and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country at the end of the product lifetime.

## Warranty

On top of the guarantee under statutory law we grant a two-year warranty on the product commencing on the date of purchase.

Within the warranty period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. Depending on availability replacement of a unit might result in a different colour or an equivalent model.

This warranty extends to every country where this appliance is supplied by Braun / Oral-B or its appointed distributor. This warranty does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The warranty becomes void if repairs are undertaken by unauthorised persons and if original Braun / Oral-B parts are not used.

This appliance is equipped with a *Bluetooth* Smart approved radio interface. Failure to establish a *Bluetooth* Smart connection to specific smart-phones is not covered under the device warranty, unless the malfunction of the *Bluetooth* interface is referable to Oral-B. The *Bluetooth* devices are warranted by their own manufacturers and not by Oral-B. Oral-B does not influence or make any recommendations to device manufacturers, and therefore Oral-B does not assume any responsibility for the number of compatible devices with our *Bluetooth* interface.

Oral-B reserves the right, without any prior notice, to make any technical modifications or changes to system implementation of device features, interface, and menu structure changes that are deemed necessary to ensure that Oral-B systems function reliably.

To obtain service within the warranty period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Customer Service Centre.

## Information

### **Bluetooth® radio module**

Although all indicated functions on the *Bluetooth* device are supported, Oral-B does not ensure 100% reliability in connection, nor feature operation consistency.

Operation performance and connection reliability are direct consequences of each individual *Bluetooth* device, software version, as well as the operating system of those *Bluetooth* devices and company security regulations implemented on the device.

Oral-B adheres to and strictly implements the *Bluetooth* standard by which *Bluetooth* devices can communicate and function in Oral-B iO toothbrushes.

However, if device manufacturers fail to implement this standard, *Bluetooth* compatibility and features will be further compromised and the user may experience functional and feature-related issues. Please note that the software on the *Bluetooth* device can significantly affect compatibility and operation.

### **Warranty replacement brush heads**

Oral-B warranty will be voided if electric rechargeable toothbrush handle damage is found to be attributed to the use of non-Oral-B replacement brush heads. Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

## Troubleshooting

Problem	Possible reason	Remedy
<b>ORAL-B APP</b>		
Toothbrush does not operate (properly) with Oral-B App.	<ol style="list-style-type: none"><li>1 Oral-B App is closed.</li><li>2 <i>Bluetooth</i> is not activated on smartphone.</li><li>3 <i>Bluetooth</i> is disabled on toothbrush.</li><li>4 <i>Bluetooth</i> connection to smartphone has been lost.</li><li>5 Your smartphone does not support <i>Bluetooth</i> 4.2 (or higher) / <i>Bluetooth</i> Smart.</li><li>6 Outdated Oral-B App.</li><li>7 Your smartphone is not compatible with the Oral-B App.</li></ol>	<ol style="list-style-type: none"><li>1 Launch Oral-B App.</li><li>2 Activate <i>Bluetooth</i> on smartphone (described in its user manual).</li><li>3 Enable <i>Bluetooth</i> via display Settings → <i>Bluetooth</i> enable.</li><li>4 Pair toothbrush and smartphone again via App settings. Keep your smartphone in close proximity when using it with toothbrush.</li><li>5 Your smartphone must support <i>Bluetooth</i> 4.2 (or higher)/<i>Bluetooth</i> Smart in order to pair toothbrush.</li><li>6 Download latest Oral-B App.</li><li>7 Check <a href="http://app.oralb.com">app.oralb.com</a> for compatible smartphones and country availability of the Oral-B App.</li></ol>

<b>TOOTHBRUSH</b>		
Unwanted toothbrush behaviour	Toothbrush was misconfigured	Perform a brush restart to restore default settings. Navigate on toothbrush display to Settings → Restart Brush (23). Press on/off button to confirm. Or use the app to perform brush restart: Devices → your brushing system → reset brush settings.
Toothbrush is not operating	<ol style="list-style-type: none"> <li>1 Battery has a very low charge level; toothbrush does not operate.</li> <li>2 Operating outside standard temperature range (message on display and light ring illumination in orange).</li> </ol>	<ol style="list-style-type: none"> <li>1 Charge at least 30 min.</li> <li>2 Always operate in standard operating temperature range (&gt; 5 °C and &lt; 35 °C). Activate toothbrush by placing it on the charger.</li> </ol>
Toothbrush operates a short stuttering sound after 2 minutes or does not operate any stuttering sound any more.	Timer has been modified/ deactivated via Oral-B App.	Use the Oral-B App to change timer settings or perform factory reset (see App Trouble Shootings).
Toothbrush does not charge.	<ol style="list-style-type: none"> <li>1 Battery is at a very low charger level (0%); it may take up to 10- 15min until display turns on.</li> <li>2 Ambient temperature for charging is out of valid temperature range (&lt; 5 °C and &gt; 35 °C).</li> <li>3 Charger may not be (properly) connected.</li> <li>4 You may be using a non-compatible / non- Oral-B charger.</li> </ol>	<ol style="list-style-type: none"> <li>1 Charge at least 30 min.</li> <li>2 Recommended ambient temperature range for charging is 5 °C to 35 °C.</li> <li>3 Ensure charger is connected to mains voltage.</li> <li>4 Ensure you are using an Oral-B charger that is compatible with your iO brush.</li> </ol>
Rattling/loud noise during brushing or turning on.	<ol style="list-style-type: none"> <li>1 Brush head may not be placed correctly onto the brush handle.</li> <li>2 Brush head may be worn out.</li> </ol>	<ol style="list-style-type: none"> <li>1 Check proper fit of brush head on handle. Small gap remains in between.</li> <li>2 Replace brush head.</li> </ol>

The rechargeable iO toothbrush contains a radio module which operates in the 2.4 to 2.48 GHz band with a maximum power of 1 mW.

Hereby, Braun GmbH declares that the radio equipment types 3786, 3796, 3795, 3758, 3758a, 3753, 3776, 3794, 3771, 3765 and 3767 are in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address:  
[www.oralb.com/ce](http://www.oralb.com/ce).

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